
PETER **BAUSCH**

Professional
Pedicure & Pedipeel

Art. 0300

Gebrauchsanweisung • Instructions for use • Mode d'emploi
Istruzioni per l'uso • Instrucciones para el uso • Gebruiksaanwijzing

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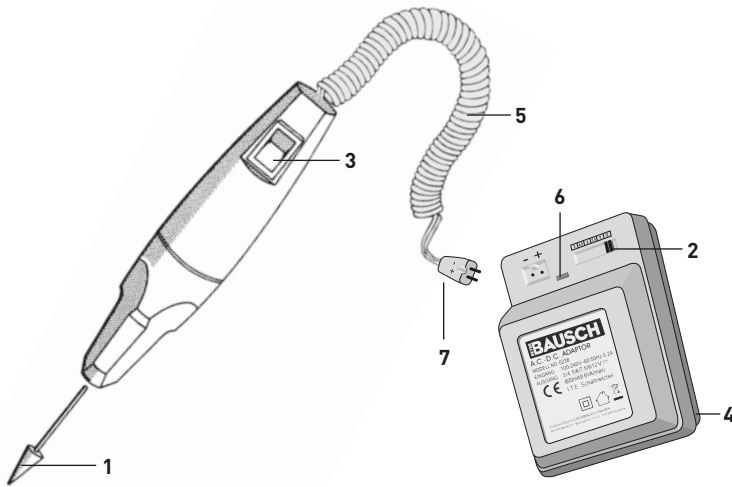1. Interchangeable tip
2. Speed regulator
3. On/Off
4. Adaptor
5. Spiral cable
6. Pilot lamp
7. Small plug

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Dear customer,
congratulations for having purchased an electric manicure and pedicure set for perfect hands and feet and thank you for having chosen our product, surely due to the quality and the design as these are our principal aims together with the satisfaction and enjoyment of our clients.

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 Please, read and follow these instructions of safety before using this appliance. Our warranty doesn't cover damages and/or injuries caused by non-observance of these safety norms and instructions. Keep this booklet of instructions in a safe place and adjoint it to the appliance if or when you should give it away.

■ IMPORTANT SAFEGUARDS

Danger! Never use the appliance in the bath, or near a water-filled hand-basin.
Risk of a fatal electrical shock!

- Do not let children use the equipment without supervision. Injury and electrical shock hazard!
- Connect the adaptor exclusively to AC mains of 230 Volts ~ 50Hz. Device may ignite spontaneously!
- For cleaning, never immerse the device in water. There's a risk of electrical shock due to penetrated water!
- Connect the appliance exclusively to the adapter supplied with it. Risk of fire!
- Never bend or squeeze the cable, nor let it get in contact with hot things. Risk of electrical shock due to damaged cable insulation!
- After each use unplug the adapter from the mains. Risk of fire due to over-heating of the transformer!
- Before each use check the appliance, in case of damage of any kind avoid using it. Injury and electrical shock hazard!
- **Attention!** Do not open or damage the device as well as the pedicure device by forcefully impacts. Injury and electrical shock hazard!

■ STARTING UP

Press the chosen tip firmly into the hole at the point of the appliance.

Connect the appliance to the adapter (4) with the small plug (7), conforming the poles (+/-). Plug in the adapter. The pilot lamp (6) indicates proper connexion to mains. Adjust the speed of rotation:

speed 1 (3 V) and 2 (4.5 V) = slow
speed 3 (6 V) and 4 (7.5 V) = middle
speed 5 (9 V) and 6 (12 V) = fast

Switch on the set at the handle!

■ USE OF THE ACCESSORIES

Hands and feet can be treated a lot easier after a warm bath (with or without additives). Nails are then more flexible, and hard skin is easier to remove.

Pointed abrasive tip: for shortening and shaping finger and toe-nails.

Cylindrical abrasive tip: to remove hard skin from fingertips, toes, heels and sole.

Cleaning and polishing brush: use this to remove skin and nail dust, or to polish the nail top. Will also clean underneath nail-ends.

Felt polishing tip: to polish the nails prior to application of lacquer. Using regularly this polisher will prevent nails to split or break.

Diamond cutter: to remove hard skin from fingertips, toes, heels and sole.

■ CLEANING AND MAINTENANCE

For cleaning just wipe the manicure/pedicure appliance with a dry, soft cloth. Don't ever put the set in water or use detergents. Except for the polishing felt and the brush, you can use water, soap and a brush to clean the other accessories. After the cleaning it is advisable to spray on some disinfectant. Let the tips dry in order to avoid rusting. The set needs no maintenance, just take good care of it.

■ REPAIRS

Caution! Electrical appliances should only be repaired by qualified technicians, as otherwise serious damage may be caused.

■ WASTE DISPOSAL

At the end of its working life the appliance must not be disposed of as urban waste. It must be taken to a special local authority differential waste collection centre for electrical and electronic products. elektronischen Geräten abgegeben werden.

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