PETER

**BAUSCH**

*Professional*  
Nail Sealer

Art. THP 1

Gebrauchsanweisung – Instructions for use – Mode d'emploi  
Istruzioni per l'uso – Instrucciones para el uso – Gebruiksaanwijzing

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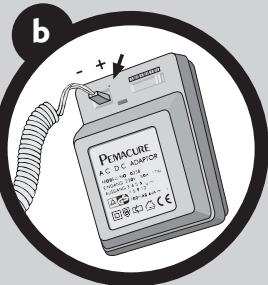
b

### GB – Nail sealer THP 1

Who'd want brittle and splitting nails?! Nails consist of 3 layers of keratin. With this patented nail sealer THP1 you can seal and reinforce them at the tip, preventing so their breaking. The heat strengthens thin and limp nails sealing together the 3 layers. It is also a very valid help in mycosis treatment.

### Starting up

Connect the sealer to the adapter (**III. b**) and plug in the adapter to the tap. When connecting the sealer to the adapter, have the poles match: plug on the adapter with the corresponding +/- poles.

The temperature is adjusted by means of the cursor:

Level 1 (3 V) and 2 (4.5 V)	= preheating of the sealing bit
Level 3 (6 V)	= sealing of very thin and limp nails
Level 4 (7.5 V) and 5 (9 V)	= sealing of brittle nails
Level 6 (12 V)	= sealing of hard nails

c

### Caution! It gets very hot! Avoid contact with the skin! Risk of burns!

**Hint:** For best performance choose the highest possible temperature you can bear and move the sealer slowly around or over the nail tip. Nails are very sensible to the heat and everyone reacts differently!

**Use:** The sealing procedure is quick and effective. When nails are long enough (**III. c**): Putting the nail tip in the groove, move the sealer slowly along the edge of the nail, back and forth 3 to 4 times. With short nails (**III. d**): Use the flat part of the sealer, moving it at the highest bearable temperature over the nail tip. Avoid contact with the skin pulling down the finger tip. The damaged nail tip is sealed and will withstand the progress of splitting and breaking. Careful! The surface of the nails is very sensible to heat. The keratin of the nails is melting gently and can no longer be damaged. The nails are protected against breakage and splitting.

**Advice:** Use sealer only on the nail tips! Don't let it come into contact with the skin, i.e. with the fingers or any other part of the body! When cutting or filing the nails, you remove the sealed part and the nail is newly exposed to brittleness and splitting. Use it as well on the toe nails (important when frequenting swimming-pools, fitness clubs, sauna, hotels, etc.).

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Peter Bausch GmbH & Co. KG